




# RevoLock®


## CLICK V2 & V3 HC Lace Kit Port A

Revised Date May 2025

RP-312A

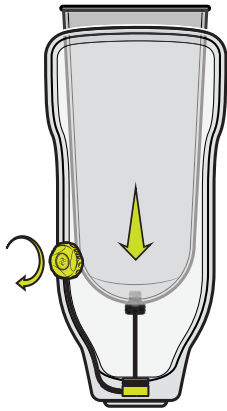

### RevoLock V2 & V3 HC Lace Kit Port A Contents:

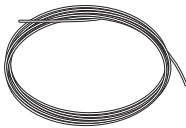

Lace



Lace Feeder

Use these instructions for devices that utilize the "A lace port" to route the lace to the distal end of the socket.



**WARRANTY:**  
www.clickmedical.co/terms/#warranty  
**U.S. Customers:**  
Contact Click Medical directly  
Help@ClickMedical.co  
Help Line: +1 970 670 7012  
**International Customers:**  
Please contact your local distributor

Click Medical products are patent protected. For full list see [www.clickmedical.co/patents](http://www.clickmedical.co/patents)

Revised 2023



### Advisory

Must be installed by a professional trained in RevoLock fabrication.  
Lace will wear. At minimum, inspect monthly and replace every 6 months.  
Loss of suspension can occur if lace fails.  
Do not use on patients with nerve issues.  
Keep lace away from open flame and sharp edges.  
Secure insert to liner before tensioning system.  
See advisories in other languages: [clickmedical.co/advisories](http://clickmedical.co/advisories)

**UK REP**

MDSS-UK RP LIMITED, 6 Wilmslow Road  
Rusholme, M14 5TP Manchester  
United Kingdom

**CH REP**

MDSS CH GmbH, Laurenzenvorstadt 61  
5000 Aarau, Switzerland

**MD**

MDSS GmbH, Schiffgraben 41  
30175 Hannover, Germany

**EC REP**




Click Medical, LLC, 1205 Hilltop Parkway, W101  
Steamboat Springs, CO 80487, USA +1-970-670-7012

For video education on Click Adjustability and RevoLock operation  
Visit [www.clickacademy.co](http://www.clickacademy.co)


### 1 Remove dial.

- Locate the metal insert.
- Insert a small flathead screwdriver between the metal insert and the dial body.
- Gently pry the reel upwards.

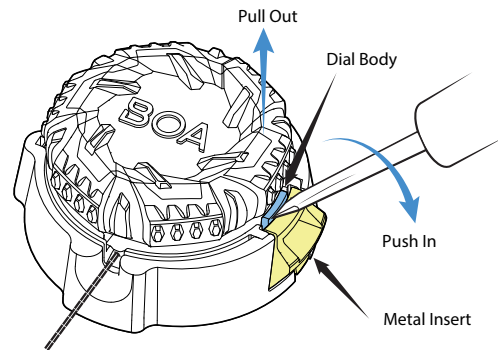


### 2 Remove the lace.


- Ensure dial is open by pulling it up.

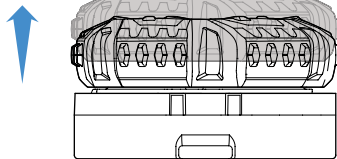

- Pull all remaining lace out of the dial.

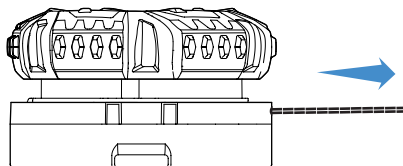


- Locate the large lace exit opening in the dial base. Hook lace loop in dial with a dental pick and pull lace out of dial.

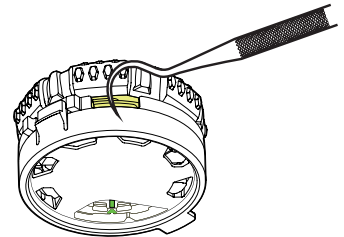


- Cut lace and pull it out of the dial. Then remove lace from the device.

Confirm that you are using the correct replacement reel by identifying the C on the bottom of the dial.

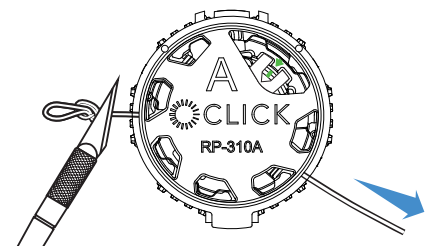


### 3 Lace the system.

Using lace feeder, route lace from the dial collar, into the tubing, through the distal housing into the socket.

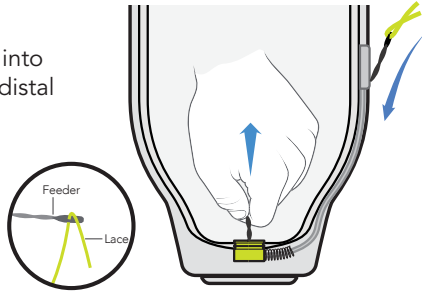


### 4 Lace the dial. a) Ready the dial for lacing.


1 Ensure dial is engaged by pushing it down.

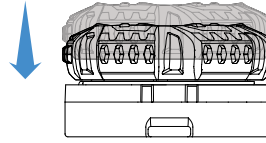


2 Turn dial over and ensure green mark is aligned with green arrow. If not, hold base and turn dial until marks align.

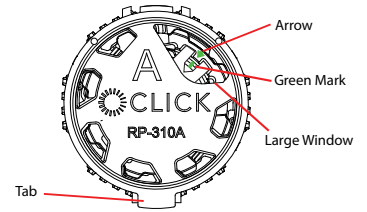

### b) Locate the lacing ports.

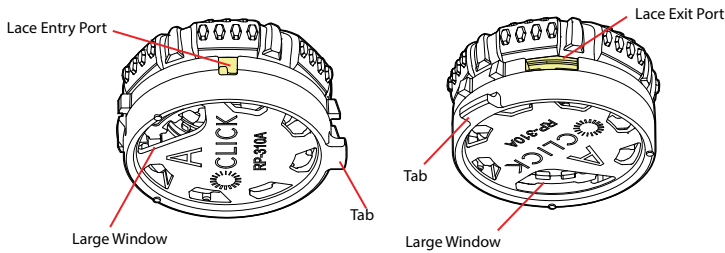


c) Feed the lace through the entry port and out the lace exit port on the opposite side of the dial.

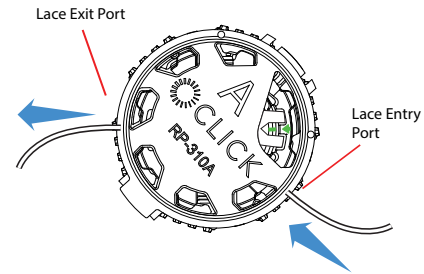

d) Tie lace with an overhand knot.


1 Make a loop.

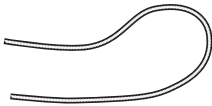

2 Tie an overhand knot with the loop.




3 Pull knot tight. Knot loop should be approximately 1cm.

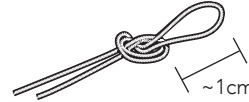

4 Trim the tail.




e) Pull loop into dial to seat knot.

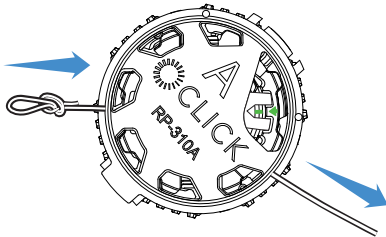

### 5 Install the dial.

a) Insert the foot of the dial into the collar.

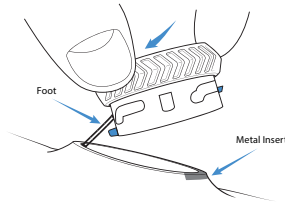


b) Press the reel firmly into the collar (you should hear a "click" when secure).




### 6 Determine lace length.

Determine how much lace is needed to extend outside of socket and attach to liner and to threaded insert.

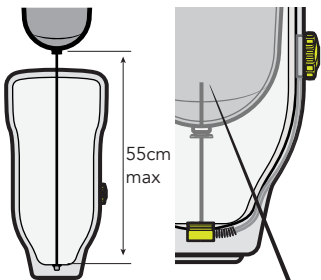

Using a lace feeder, insert lace and tie double overhand knot.

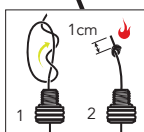


Note: Max lace length cannot exceed 55cm.

### 7 Test for function.

- Attach threaded insert to liner.
- Wind up lace.
- Pull dial out to release.
- Pull liner and lace out of socket.
- Wind up lace again.

Notes:

- Inspect lace routinely for signs of fraying or fatigue.
- Replace lace every 6 months.
- Keep lace away from open flame.

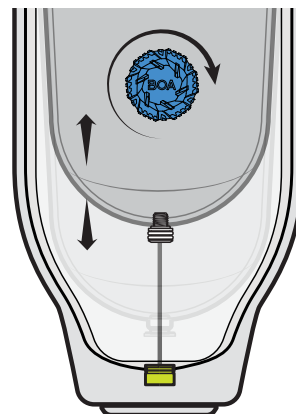

✓ Congratulations! You are finished!

